


## DVMax - Automatically Download Results from Zoasis

(Excerpts from DVMax Manual Guide 7011)

The Zoasis Lab Result function allows DVMax to automatically download your Antech Lab Results from Zoasis.com, then read downloaded Antech lab results from a specified folder and append the results to the appropriate patient record. This is available in versions 7.0.10 and higher of DVMax.

To enable the lab results reader:


- Select **Preferences** from the **File** menu.
- The **Preferences** window will appear.
- Select the **Lab Results Setup** tab.
- In the Zoasis area, select **Check for Results**.

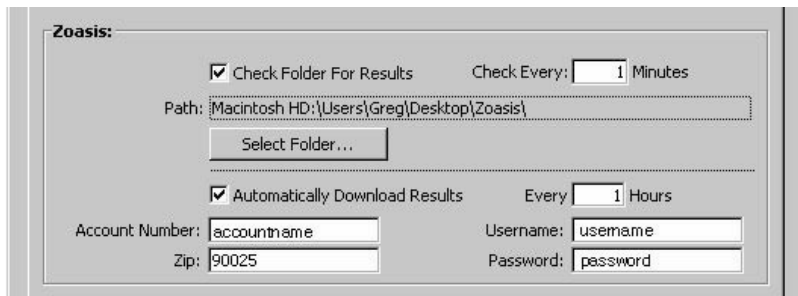

- Type a number of minutes in the **Check Every** box. This tells DVMax to check the Results folder at the specified interval.
- Type the path to the Zoasis results folder in the **Path** box, or click **Select Folder** and select the Zoasis results folder from the directory. DVMax will look for downloaded Antech lab results in this folder and place the lab results in the appropriate patient record(s).

If you would like the results to automatically download select the **Automatically Download Results** check box.

- Enter your **Zoasis ID Number, Username, Zip** and **Password**
- Click **OK** to enable checking of Zoasis lab results.

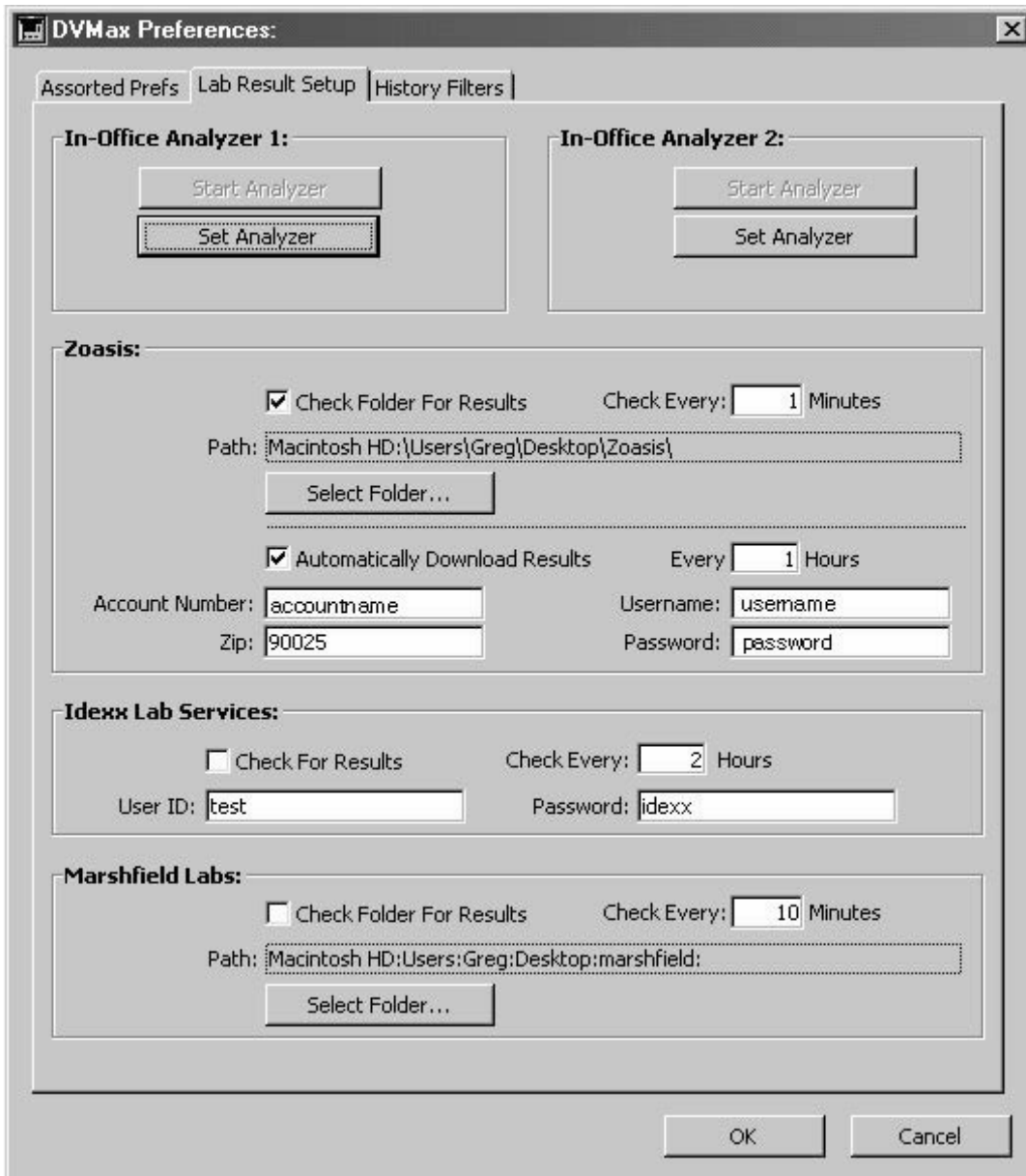

## Disable Checking for Zoasis Lab Results

To disable checking for lab results from the Zoasis folder, deselect **Check for Results** box from the Zoasis area of the Lab Result Setup tab.

